

Disposable Smoke Evacuation Forceps

Precision Grasping With Smoke Evacuation



2 - In - 1 Solution
for electrosurgery & smoke evacuation
reduces surgical plume

Hand Switch

360° Swivel

Monopolar

Bipolar

Smoke Plume Capture

Available In Monopolar
Non-Evacuation version
with Non-coated Tip



Accessories

▪ Tips Size



1.20mm
Serrated



1.20mm
Smooth



0.75mm
Serrated



0.75mm
Smooth



0.5mm
Smooth



0.1mm
Smooth

▪ Compatible Accessories

Cable and Tubing



- Cable and tubing available in 3 meters and 4.5
- Tubing and connector available in different sizes

Plug




Monopolar



Bipolar

- Compatible with all the generator such as ValleyLab, Bovie, and Conmed

Features

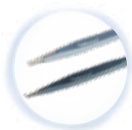

Available in Both Bipolar and Monopolar Configurations

Both Bipolar and Monopolar configurations are available catering to different needs.



Smoke Evacuation

Extracting hazardous surgical smoke at the tip of the forceps with incorporated smoke tube.

Non-Evacuation monopolar forceps is also available.



Hand Activation Control

Completely hand-activated control by squeezing and holding the buttons on the sides; Footswitch is not necessary during surgery.



Flexible Tubing and Lightweight

Reducing the operator fatigue while maximizing the flexibility by the integration of smoke tubing with 360° swivel connection.



Tips with Optional PTFE Coating

Reducing the risk of building up eschar on the tips and delivering reliable electrosurgical performance.

PTFE coating is available in Monopolar Forceps to eliminate tissue adhesion further.

Disposable Smoke Evacuation Forceps



PRODUCT CODE	PRODUCT DESCRIPTION	TIP TYPE	TIP SIZE	CABLE LENGTH (M)	PLUG TYPE
FBS001	Single Use Hand Switch Bipolar Forceps with Smoke Evacuation	Non-coated	0.1MM Smooth Tip	3M	Bipolar
FBS005	Single Use Hand Switch Bipolar Forceps with Smoke Evacuation	Non-coated	0.5MM Smooth Tip	3M	Bipolar
FBS006	Single Use Hand Switch Bipolar Forceps with Smoke Evacuation	Non-coated	0.75MM Smooth Tip	3M	Bipolar
FBS007	Single Use Hand Switch Bipolar Forceps with Smoke Evacuation	Non-coated	1.2MM Smooth Tip	3M	Bipolar
FMS001	Single Use Hand Switch Monopolar Forceps with Smoke Evacuation	Teflon coated	0.1MM Smooth Tip	3M	Monopolar
FMS005	Single Use Hand Switch Monopolar Forceps with Smoke Evacuation	Teflon coated	0.5MM Smooth Tip	3M	Monopolar
FMS006	Single Use Hand Switch Monopolar Forceps with Smoke Evacuation	Teflon coated	0.75MM Smooth Tip	3M	Monopolar
FMS007	Single Use Hand Switch Monopolar Forceps with Smoke Evacuation	Teflon coated	1.2MM Smooth Tip	3M	Monopolar