

# Telescopic Smoke Evacuation Pencil

Slide-To-Extend Technology For Complete Surgical Control



## SmokeFin

Smart · Clear · Ergonomic · Disposable

- **Enhanced Ergonomics**  
Lightweight with 360° swivel and flexible tubing.

- **Slim, Oval Triangular Grip**  
Ensures improved access and visibility.

- **Rocker Switch with Sliding Button**  
Telescopic design that slides along the housing for quick adjustments.


- **Efficiency Smoke Clearance**  
Rapid removal of surgical smoke to maintain a clear operating field.

- **Slide-to-Extend up to 60mm**  
By pushing the button forward, the housing can extend up to 60mm to reach deeper into surgical cavities.


- **Adjustable Length**  
The distance between the sliding button and electrode tip smoothly extends from 68MM to 198mm.

68MM-198MM

Hassle-free  
work-length adjust

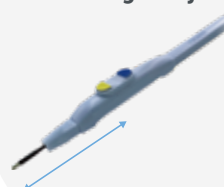

Full extension

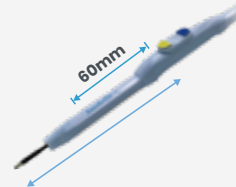

Pushing the button forward extends the electrode tip, while pulling it backward shortens the reach for precise control.

## Adjustable Reach, Efficient Smoke Evacuation

The **SmokeFin Telescopic Smoke Evacuation Pencil** features a sliding button that enables surgeons to instantly adjust the working length from 68mm to 198mm without changing the electrode. Its innovative telescopic mechanism ensures seamless adaptation during procedures, providing superior visibility, enhanced ergonomics, and effortless control.



## Accessories

### Nose Cone



- 3 sizes of nose cones fit for 63.5/69/120MM electrodes
- Transparent, equipped with radiopaque identifier (xray-detectable)

### Electrodes




- Standard 50MM  
Optional 63.5/69/120MM
- Modified blade electrode (shrink wrap insulation), coated or non-coated

### Cable and Tubing



- Standard 3 meter or 4.5 meter
- 22MM filter port adaptor
- Y-connection

### LLETZ/ LEEP electrode

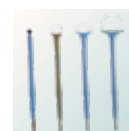

- Compatible with a range of electrode types to suit different surgical procedures

## Smoke Evacuation Pencil

- Sterile Packing
- Packing: 5pcs/Box



PRODUCT CODE	HOUSING LENGTH	ELECTRODE LENGTH	BLADE TYPE	CABLE LENGTH(M)	PLUG TYPE
PS300	195MM	50MM	Non-Coated	3M	3 Prong Plug
PS300T	195MM	50MM	Non-Sticked	3M	3 Prong Plug
PS400	195MM	63.5MM	Non-Coated	3M	3 Prong Plug
PS400T	195MM	63.5MM	Non-Sticked	3M	3 Prong Plug
PS500	195MM	69MM	Non-Coated	3M	3 Prong Plug
PS500T	195MM	69MM	Non-Sticked	3M	3 Prong Plug
PS600	195MM	120MM	Non-Coated	3M	3 Prong Plug
PS600T	195MM	120MM	Non-Sticked	3M	3 Prong Plug